CPep

Immuno Enzymometric Assay

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Attention

For North and South American Customers: Please refer to the AIA-AAM Docs on CD for the appropriate information.

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ST AIA-PACK C-Peptide

For Quantitative Measurement of C-peptide in Serum, Heparinized Plasma, EDTA Plasma or Urine

NAME AND INTENDED USE

ST AIA-PACK C-Peptide is designed for IN VITRO DIAGNOSTIC USE ONLY for the quantitative measurement of C-peptide in human serum, heparinized plasma, EDTA plasma, or urine on TOSOH AIA System Analyzers.

SUMMARY AND EXPLANATION OF TEST

C-peptide, a polypeptide 31 amino acids in length, originates in pancreatic B-cells as a metabolically inert by-product in the synthesis of insulin from proinsulin (1). Insulin and C-peptide are released from proinsulin in equimolar concentrations into the portal circulation (2). Therefore, C-peptide levels can serve as an index to insulin secretion (3,4). Where insulin secretion is diminished, as in insulin-dependent diabetes, low C-peptide levels are to be expected. Elevated C-peptide levels may result from increased B-cell activity associated with insulinomas (4-6).

Anti-insulin antibodies are commonly found in patients who have undergone insulin therapy. These circulating anti-insulin antibodies may interfere with certain immunoassays for insulin, making it difficult to use as a measure of residual B-cell activity. C-peptide measurement is, therefore, used as an alternative measurement index in this context (7). C-peptide measurement is also used as an additional means of evaluating glucose tolerance tests (8).

PRINCIPLE OF THE ASSAY

The ST AIA-PACK C-Peptide is a two-site immunoassay which is performed entirely in the ST AIA-PACK C-Peptide test cups. C-peptide present in the test sample is bound with antibody immobilized on a magnetic solid phase and enzyme-labeled monoclonal antibody in the test cups. The magnetic beads are washed to remove unbound enzyme-labeled monoclonal antibody and are then incubated with a fluorogenic substrate, 4-methylumbelliferyl phosphate (4MUP). The amount of enzyme-labeled monoclonal antibody that binds to the beads is directly proportional to the C-peptide concentration in the test sample. A standard curve is constructed, and unknown sample concentrations are calculated using this curve.

MATERIAL PROVIDED (ST AIA-PACK C-Peptide, Cat. No. 0025284)

Plastic test cups containing lyophilized twelve magnetic beads with anti-C-peptide rabbit polyclonal antibody and 100 µL of anti-C-peptide mouse monoclonal antibody conjugated to bovine alkaline phosphatase with sodium azide as a preservative.

MATERIALS REQUIRED BUT NOT PROVIDED

The following materials are required to perform C-peptide analysis using the ST AIA-PACK C-Peptide (Cat. No. 0025284) on the TOSOH AIA System Analyzers. They are available separately from TOSOH.

| Materials | Cat. No. |
|---|---------------------|
| AIA Nex-1A or AIA-21 | 0018539 |
| AIA Nex-1A or AIA-21 LA | 0018540 |
| AIA-1800 ST | 0019836 |
| AIA-1800 LA | 0019837 |
| AIA-2000 ST | 0022100 |
| AIA-2000 LA | 0022101 |
| AIA-600 II | 0019014 |
| AIA-600 II BCR | 0019328 |
| AIA-900 | 0022930 |
| AIA-360 | 0019945 |
| AIA-PACK SUBSTRATE SET II | 0020968 |
| AIA-PACK SUBSTRATE REAGENT II | |
| AIA-PACK SUBSTRATE RECONSTITUENT II | |
| AIA-PACK C-Peptide CALIBRATOR SET | 0020384 |
| AIA-PACK C-Peptide CALIBRATOR (1) | 0 ng/mL |
| AIA-PACK C-Peptide CALIBRATOR (2) | 0.5 ng/mL (approx.) |
| AIA-PACK C-Peptide CALIBRATOR (3) | 2 ng/mL (approx.) |
| AIA-PACK C-Peptide CALIBRATOR (4) | 6 ng/mL (approx.) |
| AIA-PACK C-Peptide CALIBRATOR (5) | 15 ng/mL (approx.) |
| AIA-PACK C-Peptide CALIBRATOR (6) | 33 ng/mL (approx.) |
| AIA-PACK C-Peptide SAMPLE DILUTING SOLUTION | 0020584 |
| AIA-PACK WASH CONCENTRATE | 0020955 |
| AIA-PACK DILUENT CONCENTRATE | 0020956 |
| SAMPLE CUPS | 0018581 |
| AIA-PACK DETECTOR STANDARDIZATION TEST CUP | 0020970 |
| AIA-PACK SAMPLE TREATMENT CUP | 0020971 |

Additional Requirements for AIA Nex-1A / AIA-21 only:

| | |
|---|---------|
| PIPETTE TIPS | 0018552 |
| PRELOADED PIPETTE TIPS | 0018583 |
| Additional Requirements for AIA-600 II, AIA-900, AIA-1800 and AIA-2000: | |
| PIPETTE TIPS | 0019215 |
| TIP RACK | 0019216 |
| PRELOADED PIPETTE TIPS | 0022103 |

Only materials obtained from TOSOH should be used. Materials obtained elsewhere should not be substituted since assay performance is characterized based strictly on TOSOH materials.

WARNINGS AND PRECAUTIONS

- The ST AIA-PACK C-Peptide is intended for in vitro diagnostic use only.
- Inspect the packaging and the exterior of the aluminum pouch for any sign of damage before use. If any damages are visible, contact your local TOSOH sales representative.
- Test cups from different lots or different assays should not be mixed within a tray.
- The ST AIA-PACK C-Peptide contains sodium azide, which may react with lead or copper plumbing to form potentially explosive metal azides. When disposing of such reagents, always flush with large volumes of water to prevent azide build-up.
- Human serum is not used in the preparation of this product; however, since human specimens will be used for samples and other quality control products in the lab may be derived from human serum, please use standard laboratory safety procedures in handling all specimens and controls.
- Do not use beyond the expiration date.

- The ST AIA-PACK C-Peptide has been designed so that the high dose "hook effect" is not a problem for the vast majority of samples. Serum, heparinized plasma or EDTA plasma samples with C-peptide concentrations between 30 and 90 ng/mL will read > 30 ng/mL. On the recommended 10-fold dilution, urine samples with C-peptide concentrations between 300 and 900 ng/mL will read > 300 ng/mL. The "hook effect" phenomenon may occur at C-peptide concentrations > 90 ng/mL in serum, heparinized plasma or EDTA plasma, and > 900 ng/mL in urine.
- For safe waste disposal, it is recommended that each laboratory complies with established laboratory procedures and local, state, and federal regulations.

STORAGE AND STABILITY

All unopened materials are stable until the expiration date on the label when stored at the specified temperature.

| Materials | Cat. No. |
|---|----------|
| 2-8°C: | |
| ST AIA-PACK C-Peptide | 0025284 |
| AIA-PACK C-Peptide CALIBRATOR SET | 0020384 |
| AIA-PACK C-Peptide SAMPLE DILUTING SOLUTION | 0020584 |
| AIA-PACK SUBSTRATE SET II | 0020968 |
| AIA-PACK WASH CONCENTRATE | 0020955 |
| AIA-PACK DILUENT CONCENTRATE | 0020956 |
| 1-30°C: | |
| AIA-PACK DETECTOR STANDARDIZATION TEST CUP | 0020970 |
| AIA-PACK SAMPLE TREATMENT CUP | 0020971 |

After opening the aluminum pouch, ST AIA-PACK C-Peptide test cups can be left on-board of the TOSOH AIA System Analyzers (18-25°C) for a maximum of 10 days (10 x 24 hours). When stored over night at 2-8°C, the test cups can be used for up to 30 days (30 cycles of 8 hours on board and 16 hours in the refrigerator). Once the aluminum pouch is opened, the test cups must be used within 30 days. AIA-PACK C-Peptide CALIBRATOR SET and AIA-PACK C-Peptide SAMPLE DILUTING SOLUTION should be used within 1 day and 7 days of opening or reconstituting, respectively, provided the vials or the bottles are kept tightly sealed and refrigerated at 2-8°C. They should be used after equilibrating to 18-25°C for about 30 minutes. The sample diluting solution can be used for up to 90 days provided that 1) it is used for manual dilution ONLY and 2) the bottles are kept tightly sealed and refrigerated immediately after use. Reconstituted substrate solution is stable for 3 days at 18-25°C or 7 days at 2-8°C. Working diluent and wash solutions are stable for 30 days at 18-25°C. Reagents should not be used if they appear cloudy or discolored.

SPECIMEN COLLECTION AND HANDLING

- Serum, heparinized plasma, EDTA plasma, or urine is required for the assay. Citrated plasma SHOULD NOT BE USED.
- When using serum, a venous blood sample is collected aseptically without additives. Store at 18-25°C until a clot has formed (usually 15-45 minutes), then centrifuge to obtain the serum specimen for assay.
- When using heparinized plasma or EDTA plasma, a venous blood sample is collected aseptically with designated additive. Centrifuge and separate plasma from the packed cells as soon as possible.
- For urine, collect a 24-hour urine without preservative and keep the sample refrigerated during the collections process. Record the total volume of the collection and use a well-mixed aliquot for analysis. Before assay, clear the sample by filtration or centrifugation.
- Specimen types should not be used interchangeably during serial monitoring of an individual patient. Measured concentrations may vary slightly between sample types in certain patients.
- Samples may be stored at 2-8°C for up to 24 hours prior to analysis. If the analysis cannot be performed within 24 hours, the sample should be stored frozen at -20°C or below for up to 60 days.
- Repeated freeze-thaw cycles should be avoided. Turbid serum samples or samples containing particulate matter should be centrifuged prior to testing. Prior to assay, bring frozen samples to 18-25°C slowly and mix gently.
- All urine specimens should be diluted by more than 10 fold with the AIA-PACK C-Peptide SAMPLE DILUTING SOLUTION prior to analysis. The dilution factor is identified in the test file as Sp.2. The AIA Nex-1A / AIA-21, AIA-600 II, AIA-900, AIA-1800 and AIA-2000 will automatically perform these dilutions and calculate the results.
- The sample required for analysis is 20 µL.

PROCEDURE

For the AIA Nex-1A / AIA-21, AIA-600 II, AIA-900, AIA-1800, AIA-2000 and AIA-360, please refer to their Operator's Manual for detailed instructions.


I. Reagent Preparation

- Substrate Solution**
Bring all reagents to 18-25°C before preparing the working reagent. Add the entire contents of the AIA-PACK SUBSTRATE RECONSTITUENT II (100 mL) to the lyophilized AIA-PACK SUBSTRATE REAGENT II and mix thoroughly to dissolve the solid material.
- Wash Solution**
Add the entire contents of the AIA-PACK WASH CONCENTRATE (100 mL) to approximately 2.0 L of CAP Class I water or the clinical laboratory reagent water (formally NCCLS Type I) defined by CLSI C3-A4 guideline, mix well, and adjust the final volume to 2.5 L.
- Diluent**
Add the entire contents of the AIA-PACK DILUENT CONCENTRATE (100 mL) to approximately 4.0 L of CAP Class I water or the clinical laboratory reagent water (formally NCCLS Type I) defined by CLSI C3-A4 guideline, mix well, and adjust the final volume to 5.0 L.

II. Calibration Procedure

- Calibration Curve**
The calibrators for use with the ST AIA-PACK C-Peptide have been standardized on WHO 1st IRP 84/510 (1986). The calibration curve for ST AIA-PACK C-Peptide is stable for up to 90 days. Calibration stability is monitored by quality control performance and is dependent on proper reagent handling and TOSOH AIA System maintenance according to the manufacturer's instructions. Recalibration may be necessary more frequently if controls are out of the established range for this assay or when certain service procedures are performed (e.g. temperature adjustment, sampling mechanism changes, maintenance of the wash probe or detector lamp adjustment or change). For further information regarding instrument operation, consult the TOSOH AIA System Operator's Manual.
A sample calibration curve from AIA-1800 follows and shows the algorithm used for calculating results. The concentration range of the calibration curve is displayed in 1/10 of assay range in urine specimens. The concentrations of urine specimens are calculated by multiplying the concentrations obtained on the calibration curve with the dilution factor. AIA Nex-1A / AIA-21, AIA-600 II, AIA-900, AIA-1800 and AIA-2000 will automatically calculate the concentrations of urine specimens using the dilution factor and report the results.

[6C Peg-047] Calibrator Lot. 2005504306080000
$$Y = (Ax + B\sqrt{x})^2 + C\sqrt{x} + D^2$$

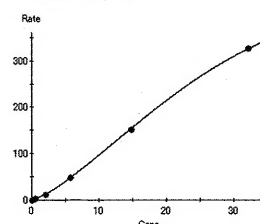B) Calibration Procedure

1. Refer to the appropriate TOSOH AIA System Operator's Manual for the procedural instructions.
 2. Verify that both the calibrator lot and concentration numbers have been correctly entered into the software.
 3. The AIA-PACK C-Peptide CALIBRATOR (1) is provided ready for use.
 4. The AIA-PACK C-Peptide CALIBRATOR (2)-(6) are lyophilized. All levels should be reconstituted with 1.0 mL of CAP Class I water or the clinical laboratory reagent water (formally NCCLS Type I) defined by CLSI C3-A4 guideline.
 5. TOSOH recommends that all calibrators be run in triplicate.
- C) Calibration Acceptability Criteria
1. The mean rate for the AIA-PACK C-Peptide CALIBRATOR (1) should be < 3.0 nmol/(L·s).
 2. Since there is a direct relationship between concentration and rate, the rate should increase as the concentration increases.
 3. The replicate values should be within a 10% range.
- D) Calibration Review and Acceptance
1. Review the calibration curve carefully, using the criteria listed above.
 2. Edit the calibration if necessary, then accept the calibration.

For further information regarding calibration, consult the TOSOH AIA System Operator's Manual.

III. Quality Control Procedure

A) Commercially Available Controls

Commercially available controls should be run at least once per day. It is recommended that at least two levels of controls, normal and abnormal, be used. Laboratory policy for this particular assay designates the following:

Control Material: _____
Frequency: _____

Lot number of control material, acceptable limits, and corrective action to be taken if controls do not meet laboratory criteria will be found in a separate quality control document maintained by the laboratory.

B) Quality Control Procedure

1. Assay quality control specimens as instructed in the specific Operator's Manual for your analyzer. In addition, refer to the TOSOH AIA System Operator's Manual for detailed instructions on defining and editing the files.
2. Quality control material to be run with this assay is defined by individual laboratory policy.

IV. Specimen Processing

A) Preparation

Serum, heparinized plasma or EDTA plasma specimens do not require dilution prior to analysis. Urine specimens should be diluted by more than 10 fold with the AIA-PACK C-Peptide SAMPLE DILUTING SOLUTION prior to analysis. Either diluted or undiluted patient specimens are placed in appropriate positions on the TOSOH AIA System Analyzers. For the AIA-360, the urine specimens must be manually diluted in sample cups and the sample cups should be placed on the AIA-360. For further information regarding sample preparation and instrument operation, refer to the SPECIMEN COLLECTION AND HANDLING section and consult the specific TOSOH AIA System Operator's Manual, respectively. Barcoded primary tubes as well as sample cups can be run on the AIA Nex-IA / AIA-21, AIA-600 II, AIA-900, AIA-1800, and AIA-2000.

B) Assay Procedure

1. Ensure a sufficient quantity of ST AIA-PACK C-Peptide test cups for the number of samples to be run.
2. Load patient samples as instructed in the Operator's Manual and proceed with analysis. Note: The AIA Nex-IA / AIA-21, AIA-600 II, AIA-900, AIA-1800 and AIA-2000 will require AIA-PACK SAMPLE TREATMENT CUPS if onboard dilutions are utilized.

PROCEDURAL NOTES

1. Lyophilized substrate must be completely dissolved.
2. Ligand assays performed by the TOSOH AIA System Analyzers require that the laboratory use water designated by the CAP as Class I or by the CLSI as the clinical laboratory reagent water. Water should be tested at least once per month and should be free of particulate matter including bacteria. The pH of the water should also be routinely tested. For further information, consult the CLSI document C3-A4, Preparation and Testing of Reagent Water in the Clinical Laboratory; Approved Guideline-Fourth Edition.
3. If a specimen C-peptide concentration is found to be greater than 30 ng/mL in serum, heparinized plasma or EDTA plasma, or 300 ng/mL in urine, the specimen should be diluted with the AIA-PACK C-Peptide SAMPLE DILUTING SOLUTION and reassayed according to the Assay Procedure. The recommended dilution for specimens containing greater than 30 ng/mL (serum, heparinized plasma or EDTA plasma) or 300 ng/mL (urine) is 1:10. It is desirable to dilute the specimen so that the diluted specimen reads between 5.0 and 15.0 ng/mL in serum, heparinized plasma or EDTA plasma, or between 50 and 150 ng/mL in urine. The dilution factor should be entered into the software. For further information on the dilution of specimens, refer to the TOSOH AIA System Operator's Manual.
4. The TOSOH AIA System Analyzers can store two different calibration curves for each analyte at one time. Therefore, up to two different lots of ST AIA-PACK C-Peptide test cups can be used during the same run.
5. If the assay specifications for this test are not ready in the system software, the specifications must be entered under test code 059.

CALCULATION OF RESULTS

The TOSOH AIA System Analyzers perform all sample and reagent handling operations automatically. The TOSOH AIA System Analyzers read the rate of fluorescence produced by the reaction and automatically convert the rate to C-peptide concentration in ng/mL.

For samples requiring dilution, the AIA Nex-IA / AIA-21, AIA-600 II, AIA-900, AIA-1800 and AIA-2000 will automatically perform dilutions and calculate results if the dilution factors are entered into the software. For detailed information regarding programming dilutions, consult the appropriate TOSOH AIA System Operator's Manual.

EVALUATION OF RESULTS

Quality Control

In order to monitor and evaluate the precision of the analytical performance, it is recommended that commercially available control samples should be assayed according to the local regulations.

The minimum recommendations for the frequency of running internal control material are:

After calibration, two levels of the internal control are run in order to accept the calibration curve.

The two levels of controls are repeated when certain service procedures are performed (e.g. temperature adjustment, sampling mechanism changes, maintenance of the wash probe or detector lamp adjustment or change).

After daily maintenance, two levels of the control should be run in order to verify the overall performance of the TOSOH AIA System Analyzers.

If one or more control sample value(s) is out of the acceptable range, it is necessary to investigate the validity of the calibration curve before reporting patient results. Standard laboratory procedures should be followed in accordance with the strict regulatory agency under which the laboratory operates.

LIMITATIONS OF THE PROCEDURE

1. For diagnostic purposes, the results obtained from this assay should be used in conjunction with other data (e.g. symptoms, results of other tests, clinical impressions, therapy, etc.).
2. Using ST AIA-PACK C-Peptide, the highest measurable concentration of C-peptide in serum, heparinized plasma or EDTA plasma specimens without dilution or in urine specimens diluted by 10 fold with the AIA-PACK C-Peptide SAMPLE DILUTING SOLUTION is 30 ng/mL, which corresponds to 300 ng/mL in original urine specimens, and the lowest measurable concentration in serum, heparinized plasma or EDTA plasma specimens without dilution or in urine specimens diluted by 10 fold with the AIA-PACK C-Peptide SAMPLE DILUTING SOLUTION is 0.2 ng/mL (assay sensitivity), which corresponds to 2.0 ng/mL in original urine specimens.

3. Although the approximate value of the highest calibrator is 33 ng/mL, the exact concentration may be slightly different. The assay specification, ASSAY RANGE HIGH, should be defined as the upper limit of the assay range, 30 ng/mL, which corresponds to 300 ng/mL in original urine specimens (dilution factor 10).
4. Although hemolysis has an insignificant effect on the assay, hemolyzed samples may indicate mistreatment of a specimen prior to assay and results should be interpreted with caution.
5. Lipemia has an insignificant effect on the assay except in the case of gross lipemia where spatial interference may occur.
6. Specimens from patients taking medicines and/or medical treatment may show erroneous results.
7. Samples containing fibrin may exhibit either falsely elevated or falsely decreased results.
8. Specimens from patients who have received preparations of mouse monoclonal antibodies for diagnosis or therapy may contain human anti-mouse antibodies (HAMA). Such specimens may show falsely elevated values when tested for C-peptide.
9. Serum concentrations typically increase post-prandially.
10. For a more complete understanding of the limitations of this procedure, please refer to the SPECIMEN COLLECTION AND HANDLING, WARNINGS AND PRECAUTIONS, STORAGE AND STABILITY, and PROCEDURAL NOTES sections in this insert sheet.

EXPECTED VALUES

Each laboratory should determine a reference interval which corresponds to the characteristics of the population being tested. As with all diagnostic procedures, clinical results must be interpreted with regard to concomitant medications administered to the patient.

Reference Ranges

The interval given here was determined in serum samples from 103 apparently healthy Asian individuals.

Serum: Reference interval = 1.1 - 3.3 ng/mL

Conversion Factors

C-peptide concentrations in this application are in units of ng/mL. Conversion to SI units of $\mu\text{g/L}$ may be made using the following equation:

$$\mu\text{g C-peptide/L} = \text{ng C-peptide/mL} \times 1.0$$

PERFORMANCE CHARACTERISTICS

ACCURACY

- a.1. Recovery (serum): Three serum pools were spiked with three different levels of C-peptide and assayed before and after spiking.

| Sample | Initial Value (ng/mL) | C-peptide Added (ng/mL) | Expected Value (ng/mL) | Measured Value (ng/mL) | Percent Recovery (%) |
|----------|-----------------------|-------------------------|------------------------|------------------------|----------------------|
| Serum A1 | 0.68 | 21.35 | 22.03 | 19.83 | 90.0 |
| | 0.68 | 10.68 | 11.35 | 10.30 | 90.7 |
| | 0.68 | 5.34 | 6.01 | 5.54 | 92.1 |
| Serum B1 | 0.21 | 21.35 | 21.56 | 19.61 | 90.9 |
| | 0.21 | 10.68 | 10.89 | 9.91 | 91.1 |
| | 0.21 | 5.34 | 5.55 | 4.93 | 88.9 |
| Serum C1 | 0.42 | 21.35 | 21.77 | 19.65 | 90.3 |
| | 0.42 | 10.68 | 11.09 | 10.10 | 91.1 |
| | 0.42 | 5.34 | 5.75 | 5.17 | 89.8 |

- a.2. Recovery (urine): One urine pool was spiked with three different levels of C-peptide and assayed before and after spiking.

| Sample | Initial Value (ng/mL) | C-peptide Added (ng/mL) | Expected Value (ng/mL) | Measured Value (ng/mL) | Percent Recovery (%) |
|----------|-----------------------|-------------------------|------------------------|------------------------|----------------------|
| Urine A1 | 4.99 | 21.35 | 26.35 | 23.59 | 89.5 |
| | 4.99 | 10.68 | 15.67 | 14.02 | 89.4 |
| | 4.99 | 5.34 | 10.33 | 9.56 | 92.5 |

- a.3. Recovery (plasma): Two plasma pools were spiked with three different levels of C-peptide and assayed before and after spiking.

| Sample | Initial Value (ng/mL) | C-peptide Added (ng/mL) | Expected Value (ng/mL) | Measured Value (ng/mL) | Percent Recovery (%) |
|-----------------------|-----------------------|-------------------------|------------------------|------------------------|----------------------|
| EDTA plasma A1 | 1.83 | 21.35 | 23.18 | 21.03 | 90.7 |
| | 1.83 | 10.68 | 12.51 | 11.71 | 93.6 |
| | 1.83 | 5.34 | 7.17 | 6.69 | 93.4 |
| Heparinized plasma A1 | 1.90 | 21.35 | 23.25 | 21.17 | 91.1 |
| | 1.90 | 10.68 | 12.57 | 12.02 | 95.6 |
| | 1.90 | 5.34 | 7.24 | 6.74 | 93.2 |

- b.1. Dilution (serum): Three serum samples containing high concentrations of C-peptide were serially diluted with the AIA-PACK C-Peptide SAMPLE DILUTING SOLUTION and assayed.

| Sample | Dilution Factor | Expected Value (ng/mL) | Measured Value (ng/mL) | Percent Recovery (%) |
|----------|-----------------|------------------------|------------------------|----------------------|
| Serum A2 | none | | 18.47 | |
| | 7.5/10 | 14.00 | 14.32 | 102.3 |
| | 5.0/10 | 9.40 | 9.79 | 104.2 |
| | 2.5/10 | 4.74 | 5.03 | 106.3 |
| | 1.0/10 | 1.91 | 2.11 | 110.4 |
| Serum B2 | none | | 20.82 | |
| | 7.5/10 | 15.83 | 16.21 | 102.4 |
| | 5.0/10 | 10.58 | 11.29 | 106.7 |
| | 2.5/10 | 5.31 | 5.81 | 109.3 |
| | 1.0/10 | 2.20 | 2.46 | 111.8 |
| Serum C2 | none | | 22.55 | |
| | 7.5/10 | 17.07 | 17.13 | 100.4 |
| | 5.0/10 | 11.46 | 11.93 | 104.1 |
| | 2.5/10 | 5.76 | 6.04 | 104.8 |
| | 1.0/10 | 2.40 | 2.58 | 107.4 |

- b.2. Dilution (urine): One urine sample containing high concentration of C-peptide was serially diluted with the AIA-PACK C-Peptide SAMPLE DILUTING SOLUTION and assayed.

| Sample | Dilution Factor | Expected Value (ng/mL) | Measured Value (ng/mL) | Percent Recovery (%) |
|----------|-----------------|------------------------|------------------------|----------------------|
| Urine A2 | none | | 18.96 | |
| | 7.5/10 | 14.24 | 14.60 | 102.5 |
| | 5.0/10 | 9.50 | 9.89 | 104.1 |
| | 2.5/10 | 4.74 | 5.15 | 108.6 |
| | 1.0/10 | 1.83 | 1.98 | 108.0 |

- b.3. Dilution (plasma): Two plasma samples containing high concentrations of C-peptide were serially diluted with the AIA-PACK C-Peptide SAMPLE DILUTING SOLUTION and assayed.

| Sample | Dilution Factor | Expected Value (ng/mL) | Measured Value (ng/mL) | Percent Recovery (%) |
|-----------------------|-----------------|------------------------|------------------------|----------------------|
| EDTA plasma A2 | none | | 22.02 | |
| | 7.5/10 | 16.61 | 17.24 | 103.8 |
| | 5.0/10 | 11.13 | 11.82 | 106.2 |
| | 2.5/10 | 5.60 | 6.05 | 108.0 |
| | 1.0/10 | 2.28 | 2.52 | 110.8 |
| Heparinized plasma A2 | none | | 21.69 | |
| | 7.5/10 | 16.37 | 16.94 | 103.5 |
| | 5.0/10 | 10.98 | 11.76 | 107.0 |
| | 2.5/10 | 5.50 | 6.01 | 109.1 |
| | 1.0/10 | 2.20 | 2.44 | 110.8 |

CORRELATION

- a. The correlation between serum (x) and heparinized plasma (y) on ST AIA-PACK C-Peptide was carried out using 50 patient specimens.

| | |
|-------------------------|------------|
| Slope | 0.978 |
| y-Intercept | -0.015 |
| Correlation Coefficient | 0.997 |
| Range of Samples | 1.0 - 11.5 |
| Number of Samples | 50 |

- b. The correlation between heparinized plasma (x) and EDTA plasma (y) on ST AIA-PACK C-Peptide was carried out using 50 patient specimens.

| | |
|-------------------------|------------|
| Slope | 0.972 |
| y-Intercept | 0.085 |
| Correlation Coefficient | 0.996 |
| Range of Samples | 1.0 - 11.4 |
| Number of Samples | 50 |

PRECISION

- a. Within run precision was determined using six controls in a total of 20 runs. Within each run, one set of duplicates per control was assayed. The mean of each duplicate was used to obtain the pooled standard deviation (SD), which was then used to calculate the coefficient of variation (CV).

| Sample | Mean (ng/mL) | Pooled SD (ng/mL) | CV (%) |
|-----------------------|-----------------|----------------------|-----------|
| Serum A3 | 4.67 | 0.0843 | 1.8 |
| Serum B3 | 8.88 | 0.1341 | 1.5 |
| Serum C3 | 18.06 | 0.2623 | 1.5 |
| Heparinized plasma A3 | 8.88 | 0.1288 | 1.5 |
| EDTA plasma A3 | 10.22 | 0.1472 | 1.4 |
| Urine A3 | 2.23 | 0.0393 | 1.8 |

- b. Total precision was determined by the duplicate assay of six controls in 20 separate runs. The means of each run were used to calculate the pooled standard deviation (SD) and coefficient of variation (CV).

| Sample | Mean (ng/mL) | Pooled SD (ng/mL) | CV (%) |
|-----------------------|-----------------|----------------------|-----------|
| Serum A3 | 4.67 | 0.1401 | 3.0 |
| Serum B3 | 8.88 | 0.2617 | 3.0 |
| Serum C3 | 18.06 | 0.5004 | 2.8 |
| Heparinized plasma A3 | 8.88 | 0.2149 | 2.4 |
| EDTA plasma A3 | 10.22 | 0.2431 | 2.4 |
| Urine A3 | 2.23 | 0.0739 | 3.3 |

SPECIFICITY

The following substances were tested for cross-reactivity. The cross-reactivity (%) is the percentage of the compound which will be identified as C-peptide. If these compounds are present in the specimen at the same concentration as C-peptide, the final result will be increased by these percentages.

| Compound | Cross-reactivity (%) |
|------------|------------------------|
| C-peptide | 100 |
| Proinsulin | 31.5 |
| Insulin | N.D. |
| | (N.D.: not detectable) |

SENSITIVITY

The minimal detectable concentration (MDC) of C-peptide is estimated to be 0.2 ng/mL, which corresponds to 2.0 ng/mL in original urine specimens (dilution factor 10). The MDC is defined as the concentration of C-peptide which corresponds to the rate of fluorescence that is two standard deviations from the mean rate of fluorescence of 20 replicate determinations of the AIA-PACK C-Peptide CALIBRATOR (1).

INTERFERENCE

Interference is defined, for the purposes of this study, with recovery outside of 10% of the known concentration of the specimen after the following substances are added to human specimens.

1. Hemoglobin (up to 410 mg/dL), free bilirubin (up to 16 mg/dL) and conjugated bilirubin (up to 17 mg/dL) do not interfere with the assay.
2. Lipemia, as indicated by triglyceride concentration (up to 1,600 mg/dL), does not interfere with the assay.
3. Ascorbic acid (up to 20 mg/dL) does not interfere with the assay.
4. Protein, as indicated by human albumin concentration (up to 5.0 g/dL), does not interfere with the assay.
5. Heparin (up to 100 U/mL) does not interfere with the assay.
6. EDTA (up to 10 mg/mL) does not interfere with the assay.

REFERENCES

1. Beischer, W., Proinsulin and C-peptide in Humans. Hormones in Normal and Abnormal Human Tissues Vol 3, pp 1-43 (1983).
2. Beyer, J., et al., C-Peptide: Its Biogenesis, Structure, Determination and Clinical Significance. Giornale Italiano di Chimica Clinica 4 (Supp.1) pp. 9-22 (1979).
3. Rubenstein, A., et al., Secretion of Proinsulin C-Peptide by Pancreatic B-cells and Its Circulation in Blood. Nature 224: 697-699 (1969).
4. Horwitz, D., et al., Circulating Serum C-Peptide. New England Journal of Medicine 295: 207-209 (1976).
5. Rubenstein, A., et al., Clinical Significance of Circulating C-Peptide in Diabetes Mellitus and Hypoglycemic Disorders. Archives of Internal Medicine 137: 625-632 (1977).
6. Rendall, M., Expanding Clinical Use of C-Peptide Radioimmunoassay. Acta Diabetologica Latina 20: 105-113 (1983).
7. Hoogwerf, B. J., Goetz, F. C., Urinary C-peptide - A Simple Measure of Integrated Insulin Production with Emphasis on Body Size, Diet, and Corticosteroids. J Clin Endocrinol Metab, 56: 60-67 (1982).
8. Kaneko, T., et al., Demonstration of C-Peptide Immunoreactivity in Various Body Fluids and Clinical Evaluation of the Determination of Urinary C-Peptide Immunoreactivity. Endocrinol. Jap., 22: 207-212 (1975).

TOSOH CORPORATION

Shiba-koen First Bldg.
3-8-2 Shiba, Minato-ku, Tokyo 105-8623
Japan
Phone :+81-3-5427-5181
Fax :+81-3-5427-5220

TOSOH EUROPE N.V.

Transportstraat 4
B-3980 Tessenderlo, Belgium
Phone :+32-13-66 88 30
Fax :+32-13-66 47 49

TOSOH BIOSCIENCE, INC.

6000 Shoreline Ct., Suite 101
South San Francisco, CA94080, USA
Phone :+1-650-615-4970
Fax :+1-650-615-0415
Phone : (800) 248-6764
Fax : (800) 685-7595

European Conformity

In vitro diagnostic medical device

Consult instructions for use


Temperature limitation

Batch code / Lot number

Manufacturer

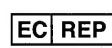Authorized representative
in the European Community

Use by date

Catalogue number
/ Part number

Supplied by

Sufficient for